

## Economy Pharmacy/Vaccine Refrigerators (+5°C)



TEPR-230-1-GD



TEPR-500-1-GD

All Thermoline Economy Pharmacy / Vaccine Refrigerators incorporate the following features:

- Internal Forced Air Flow
- Digital Temperature Display with Alarms
- Max/Min Temperature Logging
- Internal Light
- Cable Port (excludes 80 litre model)
- Door Key lock
- Frost Free Operation
- CFC Free Refrigerant
- Adjustable/Removable Shelves

### FEATURES

Three economy models of specialty pharmacy / vaccine refrigerators are available up to 500ltrs to suit situations from retail pharmacies and medical centres to universities and veterinary hospitals where pharmacy / vaccine products are required to be stored safely. A further four models are on offer up to 1275ltr capacity for increased storage requirements. (See page R-4)

Temperature is precisely controlled electronically at +5°C in ambients up to 32°C. The digital display incorporates hi (+8°C) and low (+2°C) temperature alarms, which emit visual and audible signals to warn the operator of an alarm condition. Connection to a building management system (BMS) is available as an optional extra.

A low temperature, protective, cut-out (LTPC) can be fitted as an optional extra to prevent products inside the refrigerator from freezing (particularly important for the storage of vaccines). This independent device will turn off the refrigeration system should the control fail in the cooling mode.

Cabinet construction is white colorbond exterior and molded plastic interior with curved corners for easy cleaning. Insulation is high efficiency injected polyurethane foam. All models are standard with glass self closing lockable doors. Glass doors enable the product to be viewed without loss of conditions inside the refrigerator.



TEPR-80-1-GD

APPROVED BY:



*The Quality Care Pharmacy Program (QCPP) is a quality assurance program aimed at raising the standards of service that pharmacies provide to the public. The program is based on business and professional standards developed by the Pharmacy Guild and other industry stakeholders.*

### OPTIONS

*LTPC (Low temperature protective cut out)  
Additional Shelving  
BMS (Voltage free contacts allow connection for remote monitoring of power loss, low and high temp. alarms)  
Phone Dialler (TPD-1) Analogue line only.  
POWER GUARD offering protection against Low Voltage, High Voltage Spikes*


Model	Litres	Shelves	External Dims (HxWxD)	Internal Dims (HxWxD)	Electrical	Packed Weight
TEPR-80-1-GD	80	3	89x48x53	60x38x32	3A	40kg
TEPR-230-1-GD	230	4	172x53x56	112x41x39	4A	70kg
TEPR-500-1-GD	500	5	207x61x77	142x51x51	5A	100kg

**Premium Pharmacy/Vaccine Refrigerators (+5°C)**



TPR-800-2-GD

\* Shown with optional TCR-115 and TS-1000

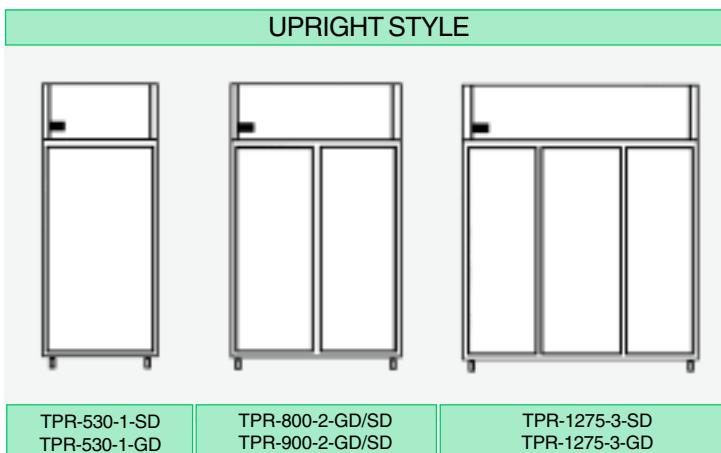

**FEATURES**

Four premium models of specialty pharmacy/vaccine refrigerators are available up to 1275 Litres to accommodate the requirements of retail pharmacies, medical centres, hospitals, universities and veterinary hospitals, where pharmacy/vaccine products are required to be stored safely. A further three economy models are on offer up to 500 Litres capacity (See Page R-3).

Temperature is precisely controlled electronically at +5°C in ambient up to 32°C. The digital display incorporates hi (+8°C) and low (+2°C) temperature alarms, which emit visual and audible signals to warn the operator of an alarm condition. Connection to a building management system (BMS) has been made easy on the premium range, as they are supplied with voltage free contacts to allow connection for remote monitoring of power loss, low and high alarms.

A low temperature, protective cut out (LTPC) is supplied as standard equipment to prevent products inside the refrigerator from freezing (particularly important for the storage of vaccines). This independent device will turn off the refrigeration system should the control fail in cooling mode.

Cabinet construction is grey colorbond exterior with interior finish of stainless steel floor and white colorbond walls. Insulation is high efficiency injected polyurethane foam. All models are available with either glass or solid self closing doors, with an optional key lock on doors for security, if required.

Premium Pharmacy / Vaccine Refrigerators incorporate the following features:

- Internal Forced Air Flow
- Digital Temperature Display with Alarms
- Max/Min Temperature Logging
- LTPC (Low Temperature Protective Cut Out)
- BMS (Voltage free contacts allow connection for remote monitoring of power loss, low and high temperature alarms)
- Internal Light
- Cable Port (on side of cabinet)
- Frost Free Operation
- CFC Free Refrigerant
- Adjustable/Removable Shelves

**OPTIONS**

- Key Lockable Doors*
- Additional Shelving*
- Phone Dialler (TPD-1) Analogue Line Only.*
- Chart Recorder (TCR-115)*
- Upgraded Alarm System 'Thermoscan 1000' (which is a battery backed up system providing a power loss and Hi Temperature. This also includes a unique 'Autotest' facility which by simply pressing a button on the control panel will warm the alarm sensor to provide a test of the alarms operation. Visual and audible signals facility provided for connection to remote monitoring system.*
- Power Guard, offering protection against Low Voltage, High Voltage Spikes.*

ORDER CODE	DIM. INT HxWxD cm	LITRES CAPACITY	DIM. EXT HxWxD cm	SHELVES	ELECTRICAL 240V 50Hz	PACKED WEIGHT Kg
TPR-530-1-GD/SD*	146x58x56	530	200x73x71	4	3A	130
TPR-800-2-GD/SD*	146x85x56	800	200x100x71	4	3A	180
TPR-900-2-GD/SD*	146x112x56	900	200x122x71	4	4A	190
TPR-1275-3-GD/SD*	146x131x56	1275	200x147x71	4	4A	210

(\*) When ordering quote **SD** for Solid Door or **GD** for Glass Door.