

Laboratory Freezers, Premium (-20°C)

AUSTRALIAN MADE



TUF-950-20-2-GD
(shown with optional TS-2000 alarm system and recorder)


* Shown with optional TCR-115 and TS-2000

UPRIGHT STYLE FAN FORCED

SD= Solid Door
GD= Glass Door



TUF-520-20-1-GD
TUF-520-20-1-SD

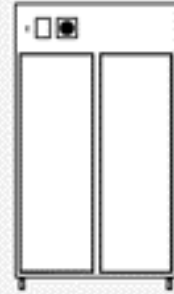

TUF-750-20-2-GD/SD
TUF-950-20-2-GD/SD

FEATURES

Thermoline have designed three -20°C premium freezers which are available in solid or double glazed glass doors. These freezers are ideal for long or short term storage of laboratory products and samples. Features include:

- Auto-Defrost feature means no more messy defrosting during critical times
- Fully adjustable open wire shelving (additional shelving available)
- Self closing front door/s
- Fan forced air distribution ensures even temperatures throughout the cabinet
- Added benefit of digital control with indication of internal conditions
- Alarms fitted for added security includes audible and visual signals with mutable function, also supplied with remote contacts at the rear for connections to BMS system (alarm probe thermally lagged to replicate product temperature and reduce false alarms)
- Easy installation - simply plug in and set temperature.

Safety Features

- ✓ Alarms fitted with rear contacts
- ✓ Over current protection
- ✓ Lockable castors

Optional Features

TCR-115, 7 day circular chart recorder

TS-2000, upgraded alarms system with built in TEST feature, battery back up and power loss notification

TPD-1, Automatic phone dialer to call when alarms are activated

POWER GUARD, offering protection against low voltage, high voltage spikes and limited lighting

ORDER CODE	DIM. INT HxWxD cm	LITRES CAPACITY	DIM. EXT HxWxD cm	SHELVES	ELECTRICAL 240V 50Hz	APPROX WEIGHT Kg.
TUF-520-20-1-GD/SD	146x58.5x56	520	200x74x71	4	4A	150
TUF-750-20-2-GD/SD	146X86X56	750	200X100X72	4	4A	270
TUF-950-20-2-GD/SD	146X109X56	950	200X147X72	4 levels	4A	380