

**Temperature and Humidity Cabinets, Refrigerated (+10°C to +60°C & 10% to 95% RH)**

AUSTRALIAN MADE

**FEATURES**

The TRH-Series of cabinets manufactured by Thermoline are a leader in the market place for accurate precise control of Temperature and Humidity. The ability to perform heat, cool, humidification and dehumidify functions makes the Thermoline TRH ideal for Pharmaceutical, Food Technology, Electronic, Bio-science and Life Science Research just to name a few. Complete in four standard capacities from 150Lt to 850Lt internal. Thermoline also has the ability to customise units to meet your application.



TRH-850  
(Shown with optional Programmable control)



Standard Controller

Eurotherm Controller



TRH-460

- Precise control of Temperatures and Humidity between +10°C and +60°C and 10% to 95% Humidity {temperature dependant}
- Digital indication of internal conditions and set values
- Microprocessor controlled with heat/cool outputs
- Internal construction of 304 stainless steel fully sealed for high humidity conditions
- Attractive PVC coated steel external surface
- Built in quiet running CFC free refrigeration system
- Fully insulated with friendly fibreglass and styrene ensures energy efficiency
- Choose between solid or glass access door/s
- Built in condensate drain system
- Fan forced air supply via long life self cooling motors delivers conditioned air throughout for precise internal conditions
- Temporal accuracy at set point better than  $\pm 0.2^{\circ}\text{C}$  & 1%Rh {unloaded chamber}
- Spatial accuracy measured to be better than  $1^{\circ}\text{C}$  & 3%Rh {unloaded chamber}
- Humidification via inbuilt vapour generating water trough with automatic water level detection and low water safety cut off

- Dehumidification via purpose built refrigeration system with air cooled condenser and electronic pulsed hot-gas by pass injection
- Heating and humidification provided through electric stainless steel heating element
- Cooling provided via evaporator coils with in air handling system and electronic switching of solenoid valves
- High grade electronic Temperature and Humidity sensor
- 50mm Monitoring or cable port access point {additional ports available}
- Open wire stainless steel fully adjustable shelves
- Convenient manual or log book holder on side of cabinet
- Easy installation with optional gravity feed water tank
- Built in condensate drain system to single exit point-allows easy cleaning of internal trough generator
- Optional ramp to set point functions
- Control error codes



TRH-150


**Temperature and Humidity Cabinets, Refrigerated (+10°C to +60°C & 10% to 95% RH)**

AUSTRALIAN MADE

- Thermoline TRH-Series Safety Devices**
- ✓ Over Current protection
  - ✓ User adjustable safety temperature cut out
  - ✓ Built in electric element safety cut out
  - ✓ Positive door closure
  - ✓ Lockable sturdy casters

- Thermoline TRH-Series Optional Accessories**
- Solid Doors or Glass doors*
  - Extended temperature range +10°C to +80°C (solid doors only) - Code: 80DEG*
  - Programmable controllers with Computer interface*
  - Retransmission of process variables*
  - Deviation alarms with remote contacts*
  - Solid viewing door over glass door (inspect product without disturbing internal conditions)*
  - Additional cable entry port holes*

**Example: TRH -150-SD-HT**  
 Temperature Range: +10.0°C to +80.0°C  
 Humidity Range: 20%RH to 95%RH

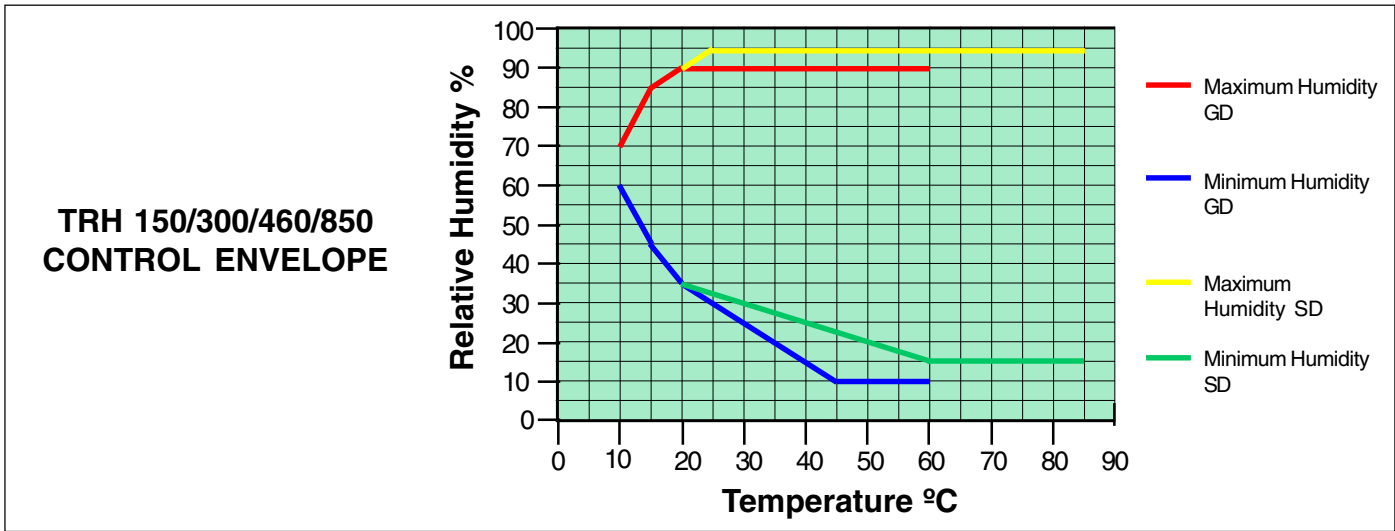

Order Code	Temperature Range °C	Humidity Range %	Internal Dims. (HxWxD) cm	Capacity Ltrs.	External Dims. (HxWxD) cm
TRH-150	10 to 60	10 to 90	64 x 50 x 50	150	135 x 66 x 76
TRH-300	10 to 60	10 to 90	99 x 60 x 50	300	175 x 76 x 76
TRH-460	10 to 60	10 to 90	114 x 65 x 65	460	190 x 81 x 91
TRH-850	10 to 60	10 to 90	114 x 117 x 65	850	190 x 133 x 91
Order Code	Nº. Shelf Positions	Nº. Shelves Supplied	Approx. Unpacked Weight Kg.	Approx. Packed Dims. (HxWxD) cm	Approx. Packed Weight Kg.
TRH-150	Multiple	2	90	157 x 87 x 91	140
TRH-300	Multiple	3	120	198 x 96 x 91	210
TRH-460	Multiple	4	150	210 x 105 x 110	230
TRH-850	Multiple	8	200	210 x 154 x 110	290
Order Code	Electrical Supply			Temperature	Humidity
TRH-150	10A 240V 50Hz	<sup>1</sup> Control Stability <sup>2</sup> Spatial Variation		±0.1 °C	±1%RH
TRH-300	15A 240V 50Hz			±1.0 °C	±3%RH
TRH-460	15A 240V 50Hz				
TRH-850	15A 240V 50Hz				

**Temperature and Humidity Cabinets, Refrigerated (+5°C to +45°C & 10% to 95% RH)**

AUSTRALIAN MADE

**FEATURES**

As well as the standard range of temperature and humidity controlled cabinets, Thermoline also has a range which includes fluorescent lighting affixed to the outside of the triple glazed door/s. These cabinets are suitable for applications where illumination of the interior is required.

The fluorescent lights are fixed to a hinged light panel which allows the light to shine through the glass door/s. This prevents any heat from the lamps affecting the conditions inside the workspace, extends the lifespan of the lamps and also allows users to easily change the lamps if required.

Conditions can also be controlled on a diurnal cycle with timers fitted as standard.



TRHL-460



Hinged door allows users to easily change the lamps if required.

**Thermoline TRHL-Series Safety Devices**


- ✓ Over Current protection
- ✓ Auto over temperature safety (Failsafe)
- ✓ Built in electric element safety cut out
- ✓ Positive door closure
- ✓ Lockable sturdy casters

**Thermoline TRHL Series Optional Accessories**

- Programmable controllers with Computer interface*
- Deviation alarms with remote contacts*
- Additional cable entry port holes*
- Extended temperature range to 60°C*
- CO2 to 3000 PPM*

**TRH 460 & 850 CONTROL ENVELOPE**

- Maximum Humidity Lights off
- Minimum Humidity Lights off
- Minimum Humidity Lights On
- Maximum Humidity Lights On

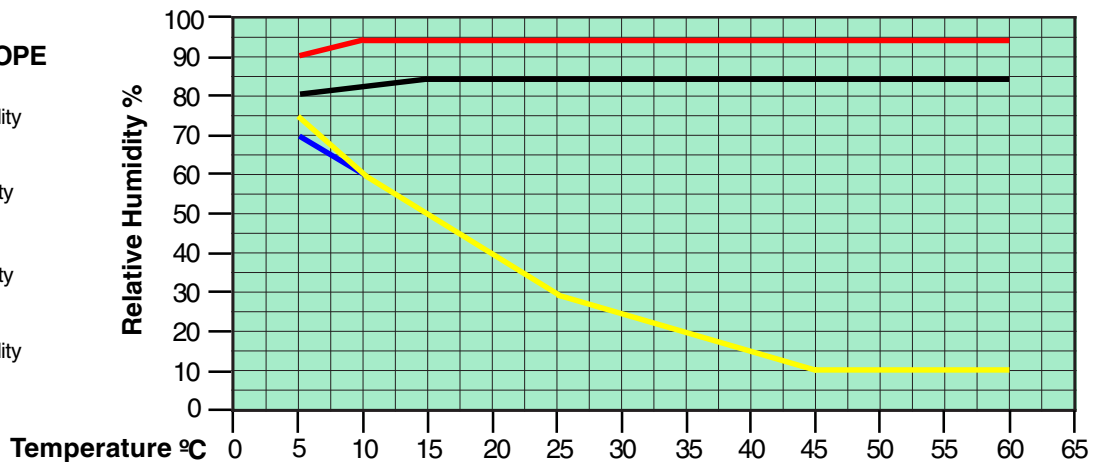

Order Code	Temperature Range °C	Humidity Range %	Internal Dims. (HxWxD) cm	Capacity Ltrs.	External Dims. (HxWxD) cm
TRHL-460	5 to 45	10 to 95	114 x 65 x 65	460	190 x 81 x 91
TRHL-850	5 to 45	10 to 95	114 x 117 x 65	850	190 x 133 x 91
Order Code	Nº. Shelf Positions	Nº. Shelves Supplied	Approx. Unpacked Weight Kg.	Approx. Packed Dims. (HxWxD) cm	Approx. Packed Weight Kg.
TRHL-460	Multiple	4	150	210 x 105 x 110	230
TRHL-850	Multiple	8	200	210 x 154 x 110	290
Order Code	Electrical Supply	Light Intensity	Temperature		Humidity
TRHL-460	15A 240V 50Hz		<sup>1</sup> Control Stability		±1%RH
TRHL-850	15A 240V 50Hz		<sup>2</sup> Spatial Variation		±3%RH

## Temperature and Humidity Cabinets, Refrigerated (+10°C to +60°C & 20% to 90% RH)

AUSTRALIAN MADE



Glass view windows

Rear perforated Air Duct

### FEATURES

Thermoline's large capacity Temperature and Humidity cabinets differ from the standard TRH-Series being designed for the long-term storage of samples where high volume capacity are needed along with precise control of all parameters.

- Two sizes 1250Lt and 2000Lt internal volume
- Precise control of Temperatures and Humidity between +10°C and +60°C and 20% to 90% Humidity {temperature dependant}
- Digital indication of internal conditions and set values
- Microprocessor controlled with heat/cool outputs
- Internal construction of 304 stainless steel fully sealed for high humidity conditions
- Attractive powder coated steel external surface
- Built in quiet running refrigeration system built below floor area




- Control error codes built in
- Fully insulated with fibreglass and styrene ensures energy efficiency
- Solid lockable access doors with internal glass viewing door
- Built in condensate drain system
- Fan forced air supply via long life self cooling motors delivers conditioned air throughout for precise internal conditions
- Temporal accuracy at set point better than  $\pm 0.2^{\circ}\text{C}$  & 1%RH (unloaded chamber)
- Spatial accuracy measured to be better than  $1^{\circ}\text{C}$  & 3%RH (unloaded chamber)
- Humidification via steam generator externally mounted under cabinet with automatic water level detection and low water safety cut off
- Dehumidification via purpose built refrigeration system with air cooled condenser and electronic pulsed hot-gas by pass injection
- Heating provided through electric stainless steel heating element and electronic pulse injection of hot gas
- Cooling provided via evaporator coils with in air handling system and electronic switching of solenoid valves
- High grade electronic Temperature and Humidity sensor
- 50mm Monitoring or cable port access point {additional ports available}
- Open wire fully adjustable stainless steel shelves
- Convenient manual or log book holder on side of cabinet
- Easy installation with optional gravity feed water tank
- Optional in ramp to set point functions

**Temperature and Humidity Cabinets, Refrigerated (+10°C to +60°C & 20% to 90% RH)**

AUSTRALIAN MADE

Order Code	Temperature Range °C	Humidity Range %	Internal Dims. (HxWxD) cm	Capacity Ltrs.	External Dims. (HxWxD) cm
TRH-1250	10 to 60	10 to 95	130 x 152 x 65	1250	202 x 170 x 82
TRH-2000	10 to 60	10 to 95	160 x 152 x 82	2000	232 x 170 x 99
Order Code	Nº. Shelf Positions	Nº. Shelves Supplied	Approx. Unpacked Weight Kg.	Approx. Packed Dims. (HxWxD) cm	Approx. Packed Weight Kg.
TRH-1250	Multiple	6	300	222 x 194 x 102	400
TRH-2000	Multiple	6	400	253 x 194 x 120	520
Order Code	Electrical Supply			Temperature	Humidity
TRH-1250	15A 240V 50Hz			<sup>1</sup> Control Stability	±0.1 °C
TRH-2000	15A 240V 50Hz			<sup>2</sup> Spatial Variation	±1.0 °C

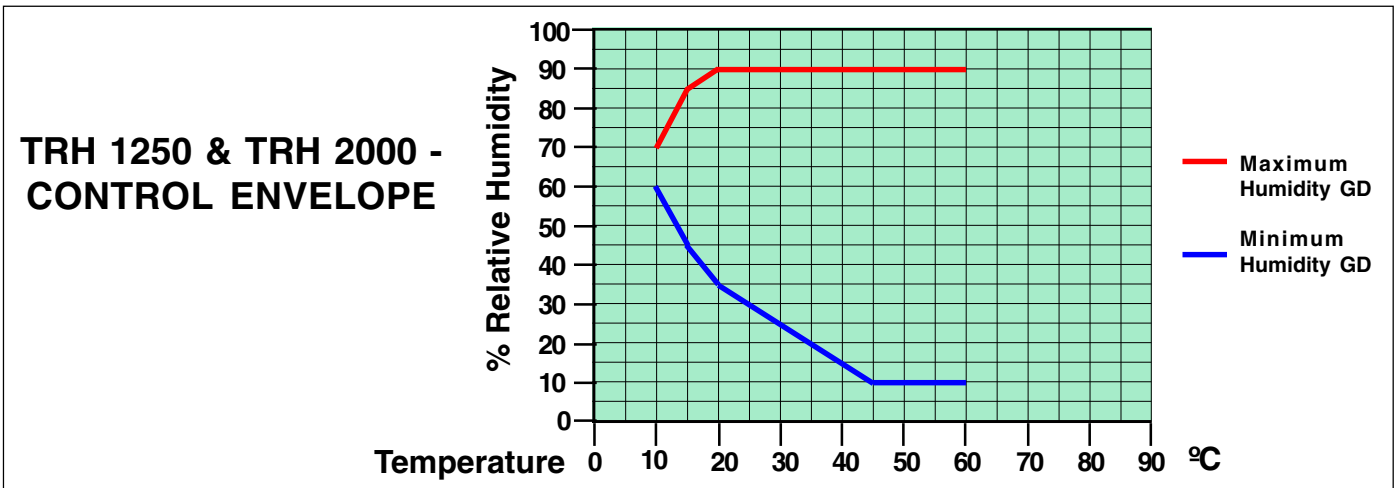

- Safety Devices**
- ✓ Over Current protection
  - ✓ Auto over temperature safety (Failsafe)
  - ✓ Built in electric element safety cut out
  - ✓ Positive door closure
  - ✓ Lockable sturdy casters

- Thermoline TRH-Series Optional Accessories**
- Programmable controllers with Computer interface*
  - Retransmission of process variables*
  - Deviation alarms with remote contacts*
  - Additional cable entry port holes*