

# SHAKING THAW BATH


(AMBIENT +5°C TO +45°C)

Image: TSBT-20 with raised rim and clear lid



## Thermoline Shaking Thaw Bath with Shaking Tray

- Adjustable temperature range up to 45°C
- Adjustable timer function (1 to 999 min)
- Accurate digital PID temperature controller
- Easy tray removal for cleaning
- Raised rim to reduce splashing
- Durable stainless steel construction
- Drain fitted standard
- Adjustable shaking stroke

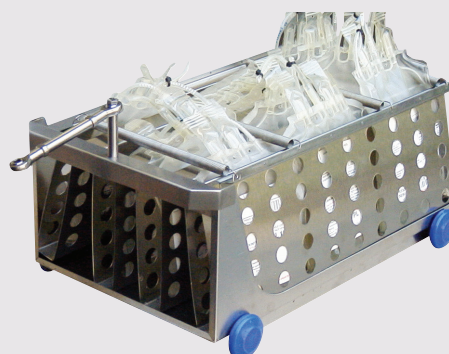

Specialty thaw tray (shown only as an example with bags)

Head Office Phone (02) 9604 3911  
Head Office Fax (02) 9725 1706

Email [info@thermoline.com.au](mailto:info@thermoline.com.au)  
Web [www.thermoline.com.au](http://www.thermoline.com.au)

 **Thermoline**  
S C I E N T I F I C  
H I G H E R S T A N D A R D S I N S C I E N C E

# PRODUCT DETAILS

This digital control shaking bath has been specially designed by Thermoline for thawing frozen products. This bath has a large capacity tray which can accommodate twelve 500ml frozen bags. The dual digital displays allow for quick and easy viewing of set temperature, actual bath temperature and allows users to easily change the operating time.

A few key features include:

- Fast Thawing time reduced to 12 minutes on the adjustable timer (full load 30 minutes)
- Unique thaw tray allows for quiet operation and perforated sides, allow the centre bags to be thawed at the same time
- Locking springs keeps all bags in the correct position

Safety Devices
Over temperature cut out set at +45°C
Auto over temperature safety (Failsafe)
Audible time out signal
Control error codes
Manual Reset Safety Device (MRST)

## SPECIFICATIONS

Model	TSBT-20
Capacity	20 litres
Tray Dimensions (LxWxH mm)	370x270x50
External Dimensions (WxDxH mm)	660x330x290
Temperature Range	Ambient +5°C to 45°C (37°C Recommended)
Temperature Stability	±0.25°C
Temperature Control	Digital PID
Weight	22kg

Thermoline Scientific have been manufacturing and distributing high quality laboratory and scientific testing equipment since 1970. Over this time, Thermoline has grown to be a leading brand in the science industry, with our equipment being used in small and large Hospitals, Universities and Research Laboratories across Australia and the Asia Pacific region.

We're proud to say that we are 100% Australian owned and operated.

**Head Office Phone** (02) 9604 3911

**Head Office Fax** (02) 9725 1706

**Email** [info@thermoline.com.au](mailto:info@thermoline.com.au)

**Web** [www.thermoline.com.au](http://www.thermoline.com.au)