
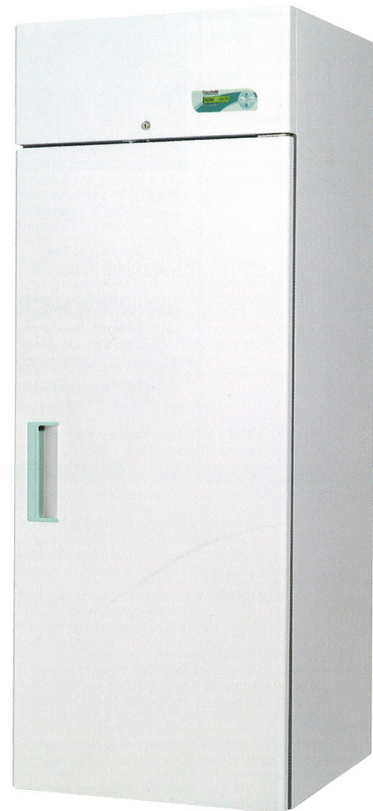


UPRIGHT SUPERARTIC FREEZERS

(-20°C TO -40°C WITH AUTO DEFROST)

Image: SUPERARTIC 250 and SUPERARTIC 600

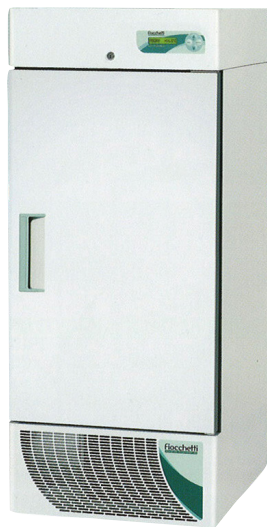

Upright Superartic Freezers

- Reversible front door to suit your location
- Lockable doors
- Adjustable shelves or optional drawers
- Forced air circulation for quick freezing
- Advanced automatic defrosting and condensate disposal
- BMS contacts to connect to a Building Management System

Head Office Phone (02) 9604 3911
Head Office Fax (02) 9725 1706

Email info@thermoline.com.au
Web www.thermoline.com.au

 **Thermoline**
S C I E N T I F I C
HIGHER STANDARDS IN SCIENCE

PRODUCT DETAILS

These two new upright freezers in the Superartic series have been specially designed for products which need to be stored at ultra low temperatures. Both freezers have an adjustable temperature range from -20°C to -40°C.

The cabinets are constructed internally with stainless steel and externally from white colorbond steel. Polyurethane foam insulation is used to help maintain the internal temperature. Heavy duty castors are used to keep the freezers stable when they are being moved. Cabinets constructed with stainless steel exterior are available at an extra cost.

Both models feature the latest controller technology with backlit alphanumeric LCD display for easy temperature adjustments and a password protect on/off option which prevents set point modifications without authorisation. Audible and visual alarms including recording of high/low temperature, power failure, door open and dirty condensor alarms notify the user should a problem occur during operation.

Safety Features

Alarms fitted standard

Over current protection

Lockable castors

Optional Accessories

Circular chart recorder

SPECIFICATIONS

Model	Superartic 250	Superartic 600
Capacity (litres)	250	600
Temperature Range	-20°C to -40°C with auto defrost	
External Dimensions (HxWxD mm)	1580x670x830	2080x720x760
Shelves	4	7
Electrical	230V 50Hz 5A	
Weight (kg)	150	280

Thermoline Scientific have been manufacturing and distributing high quality laboratory and scientific testing equipment since 1970. Over this time, Thermoline has grown to be a leading brand in the science industry, with our equipment being used in small and large Hospitals, Universities and Research Laboratories across Australia and the Asia Pacific region.

We're proud to say that we are 100% Australian owned and operated.

Head Office Phone (02) 9604 3911

Head Office Fax (02) 9725 1706

Email info@thermoline.com.au

Web www.thermoline.com.au