


REFRIGERATOR & FREEZER COMBO

(DOUBLE DOOR MODELS)

Image: LABOR 4002PC and VISION 400 2T

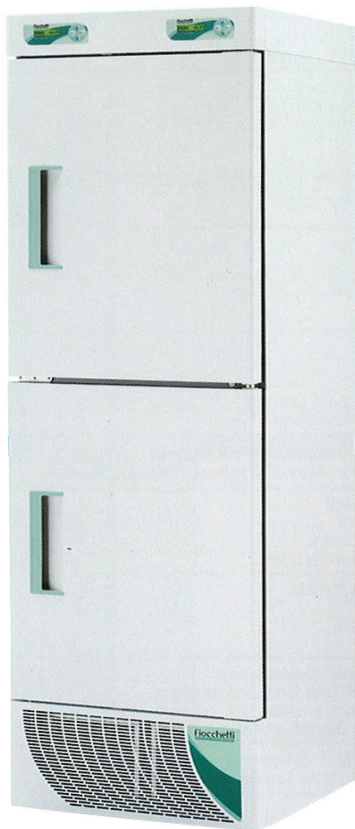

Two Temperature Refrigerators and Freezers

- Optional slide out aluminium drawers with height adjustment
- Adjustable slide out shelves (all models)
- Bottom solid door only and top glass or solid door (Labor models)
- Top and bottom double glazed glass doors (Vision model only)
- Lockable self closing door/s (all models)
- Internal lights (all models)
- Automatic condensate disposal (all models)
- Lockable castors (all models)
- BMS contacts to connect to a Building Management System (all models)

Head Office Phone (02) 9604 3911

Head Office Fax (02) 9725 1706

Email info@thermoline.com.au

Web www.thermoline.com.au

PRODUCT DETAILS

These two temperature refrigerator freezer cabinets are ideal for work areas with limited space available. Both Labor and Vision models have two independently controlled chambers - top section controllable from 0 to +15 degrees and the bottom from -15 to -20 degrees. Both sections have front opening doors which are available in both solid and glass for the Vision model, and top only glass or solid in the Labor model.

Both Labor and Vision models feature the latest controller technology with backlit alphanumeric LCD display and many unique functions including an advanced defrosting feature that allows efficient defrosting cycles which reduces temperature variations and power consumption.

These cabinets are constructed internally and externally from white colorbond steel. All models have internal rounded back corners and a dished bottom. Heavy duty castors are used to keep the cabinets stable when they are being moved.

Vision Model - TGA Approved

The Vision 400 2T model is approved by the TGA to be used as a medical device in hospitals and medical centres throughout Australia. This model has two 200 litre sections each with double and triple glazed doors. The DMLP monitor and printer is supplied standard for the top refrigerator section only. Both top and bottom sections have separate BMS monitoring.

Door Configurations
LABOR Model : Two solid doors, or bottom solid door and glass top door
VISION Model : Top and bottom glass doors

SPECIFICATIONS

Model	LABOR 2802PC	LABOR 4002PC	VISION 4002T
External Dimensions (WxDxH mm)	520x665x1650	600x630x1950	600x655x1950
Temperature Range (Top)	+2°C to +15°C		
Temperature Range (Bottom)	-15°C to -20°C		
Capacity (Top) (litres)	140	200	
Capacity (Bottom) (litres)	140	200	
Doors	Upper and lower front opening doors		
Shelves (Top and Bottom)	3 Slide out shelves		
Key Locks	Yes		
Wheel Kit	Yes		
Evaporation Tray	Yes		
DMLP Controller	No		Top section only
DMLP Printer	No		Top section only
BMS Connection	Yes		
Weight (kg)	160	180	180

Thermoline Scientific have been manufacturing and distributing high quality laboratory and scientific testing equipment since 1970. Over this time, Thermoline has grown to be a leading brand in the science industry, with our equipment being used in small and large Hospitals, Universities and Research Laboratories across Australia and the Asia Pacific region.

We're proud to say that we are 100% Australian owned and operated.

Head Office Phone (02) 9604 3911
Head Office Fax (02) 9725 1706

Email info@thermoline.com.au
Web www.thermoline.com.au

